IMPLASMILE

Impla Smile Sp. z o.o. Sp. k.

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03-885 Warszawa

tel: 22 299 21 99

kom: 790-502-199

790-552-199

e-mail: biuro@implasmile.pl

www.implasmile.pl

KISPLANT

KUWOTECH IMPLANT SYSTEM

Geo
GeoMedi

New York
KERATOR
Overdenture Attachment Systems

Medibrex

Ti-OSS[®]

KUWOTECH
IMPLANT SYSTEM

Welcome at KUWOTECH

Kuwotech was established in 2007 with the purpose of leading the advanced medical industry and has manufactured key dental materials and components based on precision process technology and started to develop, produce and sell Dental Implant and zirconia ceramic biomaterial for dental use through continuous industrial-educational joint research development. Kuwotech does our best for creating value in the future of human society with the excellent products and services based on taking care of each other, persons of talent and technology. Kuwotech, that is an enterprising company based on the spirit of perpetual challenges and future technology skills, promises to ensure your valuable future with even when facing upcoming new challenges and changes and we will be a partner who delivers happiness and help customers attain their dreams

CEO / Jeong Cheol woong

A handwritten signature in black ink, appearing to be 'Jeong Cheol woong'.

History

- 2007**
 - Establishment of MegaTech Co., Ltd.
 - Relocation of Headquarter
 - MOU on the investment with Gwangju Metropolitan City
 - Approved as medical equipment manufacturing company (No. 2603))
- 2008**
 - Establishment of R&D Center (No. 2008210228)
- 2009**
 - Change the company name (Kuwotech Co., Ltd.)
 - Quality management system ISO 9001:2000 Certified (QMS-1871)
 - Capital Increase (Capital : 1,887,500,000 KRW)
 - Ministry of knowledge economy, parts & materials professional company certified (No. 6962)
- 2011**
 - Venture business certified (No. 20110102162)
 - Ministry of knowledge economy, Regional star business awarded
 - Technology innovation small business (INNO-Biz) certified
 - Completed new factory & relocation of headquarter (50 Chemdan venture so-ro 37beon-gil, Buk-gu, Gwangju, 61003, Korea)
- 2012**
 - Headquarter of Gwangju women new job creation, selected as a women-friendly company
 - Citation as the job creation excellent company awarded, Gwangju metropolitan city (No.471)
 - Employment excellent company certified (No. 2012-8)
- 2013**
 - ISO 13485:2012 certified (Zirconia Ceramic Block)
 - CE certified (Zirconia Ceramic Block, 2195-MED-1302002)
 - FDA certified (Zirconia Ceramic Block)
 - IP star business certified (No. 2013-07)
 - Citation as the contribution to small business development awarded
- 2014**
 - CE certified (T-Cleana)
 - Citation minister of ministry of trade, industry and energy
- 2015**
 - KFDA certified (Dental Implant - 'KISPLANT')
 - Designated as a high technology company
- 2016**
 - ISO 13485:2012 certified (Dental Implant - 'KISPLANT')
 - CE certified (Dental Implant - 'KISPLANT')
- 2017**
 - Development of TREATON minocycline ointment


KISPLANT

KUWOTECH
IMPLANT
SYSTEM

Platform Switching

Conical Sealing & Internal Hex Connection

Allows for biologic width, conducting soft and hard tissue growth and reducing bone loss at osseointegrative surfaces. 11° Morse taper is self-locking to enhance the connection between fixture and abutment. Anti-rotational, secure connection between fixture and abutment.

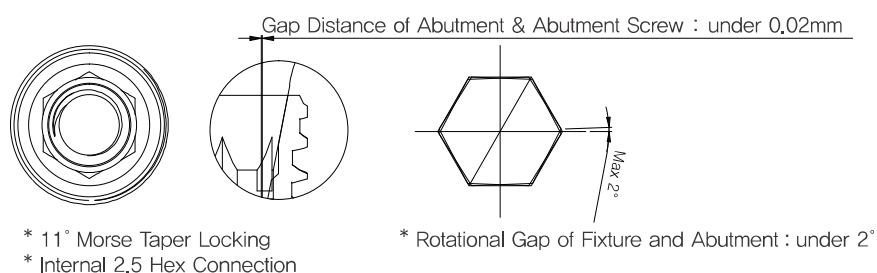
Taper-straight Design for Stable placement

Unique Macro Thread

- The design enables progressive bone condensing, gentle ridge expansion, and minimized shear force production.
- Increases the surface area in contact with bone increasing initial stability

Apex with Cutting Edge

- Enhanced self tapping based on cutting Edge
- Easy start with sharp Apex
- Reduces risk of damage to the bone

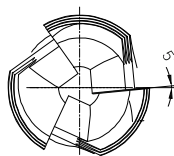Excellent soft tissue response

S-Line Abutment

Shoulder Margin for Esthetics & Zirconia Crown

Higher Compressive Strength


* Wall Thickness

	A사 Ø3.6	B사 Ø3.6	KISPLANT(Ø3.6)			A사 Ø4.0	B사 Ø4.0	KISPLANT(Ø4.0)	
			BT	MT				BT	MT
A	0.201	0.341	0.320	0.320	A	0.296	0.476	0.402	0.422
B	0.056	0.197	0.201	0.252	B	0.173	0.321	0.331	0.430
C	0.248	0.324	0.521	0.530	C	0.369	0.466	0.590	0.600

Fast & Safe osseointegration

Special S.L.A. Surface

The surface roughness of Kuwotech Co., Ltd. are achieved by sandblasting with Alumina, followed by etching with a combination of acids, a process which also creates a greater envelope surface and improves osseointegration.

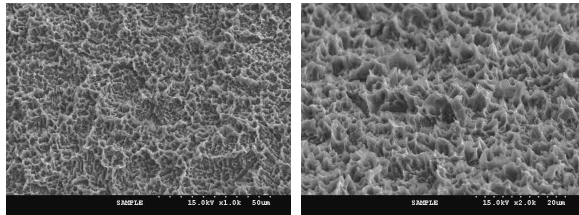Auto Etching Machine


All products of Kuwotech Co., Ltd. are produced by Auto Etching Machine which Products uniform surface treatment products.

Quality Control

All products of Kuwotech Co., Ltd. are subject to strict quality control in all process of production.

Fixture

BT(Body Thread)

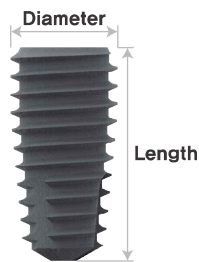- Body Thread Type for reducing load during placing
- Taper-Straight Design for stable placement
- Easy entry due to Sharp Apex
- Enhanced Self Tapping function by cutting edge
- Enclosed cover screw

Ø 3.5

Length	7.0	8.5	10.0	11.5	13.0	15.0
Ref.C	BT3507	BT3508	BT3510	BT3511	BT3513	BT3515

Ø 4.0

Length	7.0	8.5	10.0	11.5	13.0	15.0
Ref.C	BT4007	BT4008	BT4010	BT4011	BT4013	BT4015

Ø 4.5

Length	7.0	8.5	10.0	11.5	13.0	15.0
Ref.C	BT4507	BT4508	BT4510	BT4511	BT4513	BT4515

Ø 5.0

Length	7.0	8.5	10.0	11.5	13.0	15.0
Ref.C	BT5007	BT5008	BT5010	BT5011	BT5013	BT5015

Cover Screw

Ref.C KCS

- Used a Hand Driver (1,2HEX)
- Recommend Torque : 5~8Ncm

Ø 5.5

Length	7.0	8.5	10.0	11.5	13.0	15.0
Ref.C	BT5507	BT5508	BT5510	BT5511	BT5513	BT5515

Ø 6.0

Length	7.0	8.5	10.0	11.5	13.5	15.0
Ref.C	BT6007	BT6008	BT6010	BT6011	BT6013	BT6015


Fixture. BT(Body Thread)

BT(Body Thread)

- Body Thread Type for reducing load during placing
- Taper–Straight Design for stable placement
- Easy entry due to Sharp Apex
- Enhanced Self Tapping function by cutting edge
- Enclosed cover screw

Cover Screw


- Used a Hand Driver (1.2HEX)
- Recommend Torque : 5~8Ncm

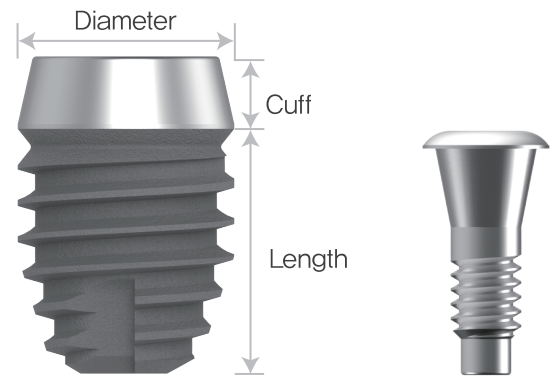Ref.C | KCS

Ø 4.5

5.0 6.0

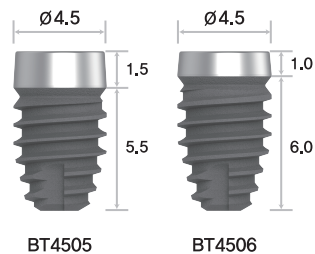
BT4505

BT4506

Ø 5.0

5.0 6.0

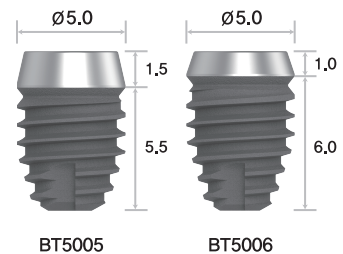
BT5005

BT5006

Ø 5.5

5.0 6.0

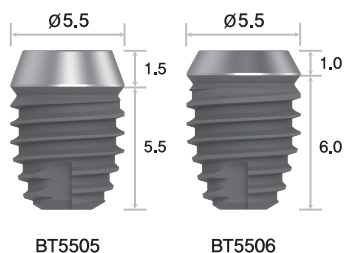
BT5505

BT5506

Ø 6.0

5.0 6.0

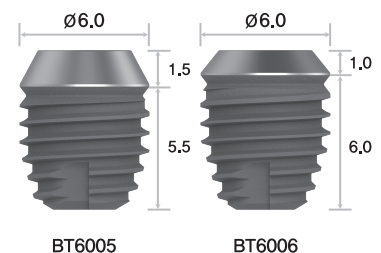
BT6005


BT6006

Drilling Sequence for Short Implant


BT4505 Final Drill : Ø4,0 Dense Drill

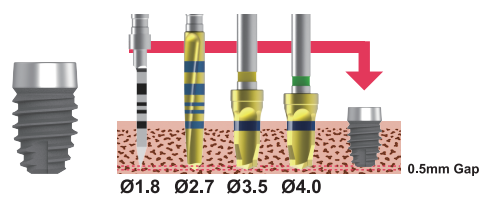
BT4506 Final Drill : Ø4,0 Dense Drill

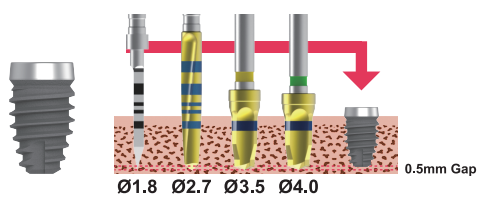
BT5005 Final Drill : Ø4,5 Dense Drill

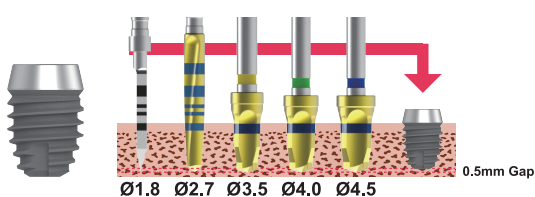
BT5006 Final Drill : Ø4,5 Dense Drill

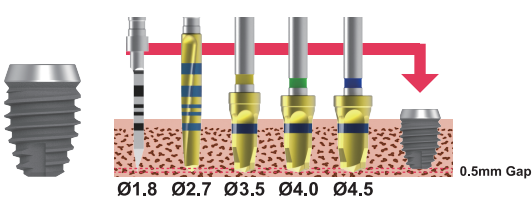
BT5505 Final Drill : Ø5,0 Dense Drill

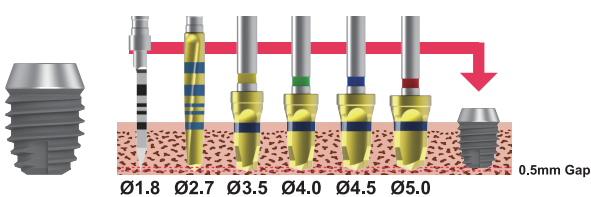
BT5506 Final Drill : Ø5,0 Dense Drill

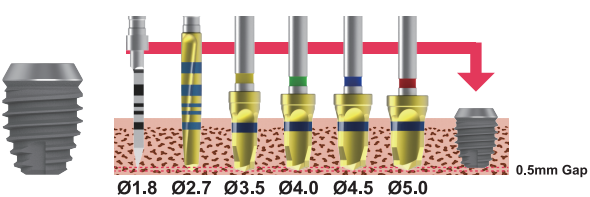
BT6005 Final Drill : Ø5,5 Dense Drill

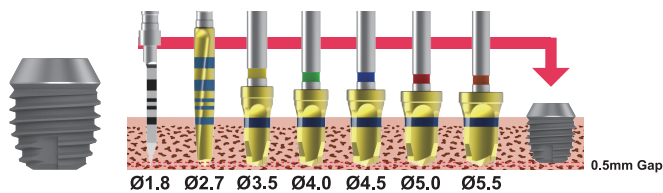BT6006 Final Drill : Ø5,5 Dense Drill

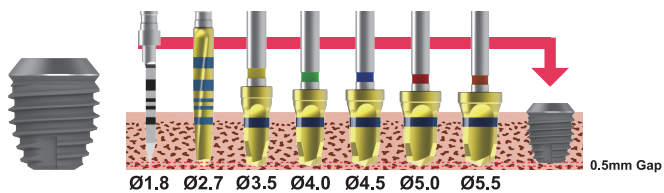MT(Micro Thread)

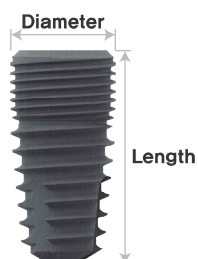- Micro double thread
- Taper–Straight Design for stable placement
- Easy entry due to Sharp Apex
- Enhanced Self Tapping function by cutting edge
- Enclosed cover screw

Ø3.5

Length	7.0	8.5	10.0	11.5	13.0
					
Ref,C	MT3507	MT3508	MT3510	MT3511	MT3513

Ø4.0

Length	7.0	8.5	10.0	11.5	13.0
					
Ref,C	MT4007	MT4008	MT4010	MT4011	MT4013

Ø4.5

Length	7.0	8.5	10.0	11.5	13.0
					
Ref,C	MT4507	MT4508	MT4510	MT4511	MT4513

Ø5.0

Length	7.0	8.5	10.0	11.5	13.0
					
Ref,C	MT5007	MT5008	MT5010	MT5011	MT5013

Ø5.5

Length	7.0	8.5	10.0
			
Ref,C	MT5507	MT5508	MT5510

Ø6.0

Length	7.0	8.5	10.0
			
Ref,C	MT6007	MT6008	MT6010

Healing Abutment

- Used a Hand Driver (1.2HEX)
- Recommend Torque : 5~8Ncm

Ø4.5

Height	3.0	4.0	5.0	6.0	7.0
Ref.C	HA4530	HA4540	HA4550	HA4560	HA4570

Ø5.5

Height	3.0	4.0	5.0	6.0	7.0
Ref.C	HA5530	HA5540	HA5550	HA5560	HA5570

Ø6.5

Height	3.0	4.0	5.0	6.0	7.0
Ref.C	HA6530	HA6540	HA6550	HA6560	HA6570

Ø7.5

Height	3.0	4.0	5.0	6.0	7.0
Ref.C	HA7530	HA7540	HA7550	HA7560	HA7570


Cover Screw

- Used a Hand Driver (1.2HEX)
- Recommend Torque : 5~8Ncm

Abutment

Fine Abutment

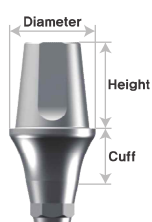- Recommend Torque : 35Ncm
- Used a Hand Driver(1.2HEX)
- Enclosed a Abutment screw

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
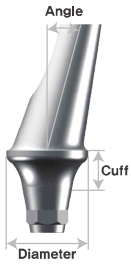
Impression Coping

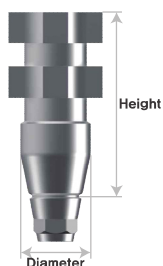- Used a Hand Driver (1,2HEX)
- Used for Closed Tray & open tray pick up (Transfer Coping) Impression Technique
- Enclosed a Guide pin

Pick up

Diameter	4.5	5.5
		
Height		
12	IPH4512	IPH5512
16	IPH4516	IPH5516

	
Guide Pin	GPP12 GPP16

Transfer

Diameter	4.5	5.5
		
Height		
12	ITH4512	ITH5512
16	ITH4516	ITH5516

	
Guide Pin	GPH12 GPH16

Lab Analog


Type	Small	Regular/Wide
		
Color	Green	Titan
Ref.C	KLA350	KLA400

- Replicates the Fixture
- Small type Lab analog for Ø3.5 Fixture
- Regular & Wide type Lab analog for all Fixture sizes except Ø3.5 Fixture

KISPLANT SYSTEM

*set (Scanbody/Multibase/CMFit Abutment + screw)

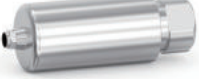
SC-Scanbody(Model)	GMS-Scanbody(Intra-Oral)	Single LL-Multibase	Bridge LL-Multibase	Single LS-Multibase
				
Ø10 CMFit Abutment	Ø14 CMFit Abutment	3D-Analog	Screw	Driver (Ncm)
				 1.25Hex 25-30Ncm

Description-Library code		0201 - DER	
Product	Connection Platform Ø	Universal	Ø3.5/4.0/4.5
SC-Scanbody			SC-DER
GMS-Scanbody(Short)			GMS-DERS
GMS-Scanbody(Long)			GMS-DERL
Ti Base	GH : 1.3		LL2-DER13-H2
	GH : 2.0		LL2-DER20-H2
	GH : 3.0		LL2-DER30-H2
	GH : 4.0		LL2-DER40-H2
Ti Base	GH : 1.3		LL2-DER13-N2
	GH : 2.0		LL2-DER20-N2
	GH : 3.0		LL2-DER30-N2
	GH : 4.0		LL2-DER40-N2
LS-Multibase / Single			LS2-DER-H2
LS-Multibase / Bridge			LS2-DER-N2
CMFit Abutment(S) / Ø10			CMFit10-DER-H2
CMFit Abutment(B) / Ø10			CMFit10-DER-N2
CMFit Abutment(S) / Ø14			CMFit14-DER-H2
CMFit Abutment(B) / Ø14			CMFit14-DER-N2
3D-Analog			3D-DER
Screw			DERC

Design concept

GMS Scanbody (Intra-oral)

- Biocompatible Titanium (Ti-6Al-4V ELI) from CARPENTER Co. guarantees durability which can be autoclavable and reusable.
- special coating allows to scan without any scan spray.
- provides with two different sizes. – Short(7.5mm) / Long(11.5mm)
- consistent performance and visible in x-rays.
- Various Implant system available.
- Highly precise Geo library is provided for CAD s/w. (3shape, exocad, Dental Wings)
- Customized Abutment Design with CNC or Premilled abutment is available.
- Zirconia prosthesis design with Geo Link.

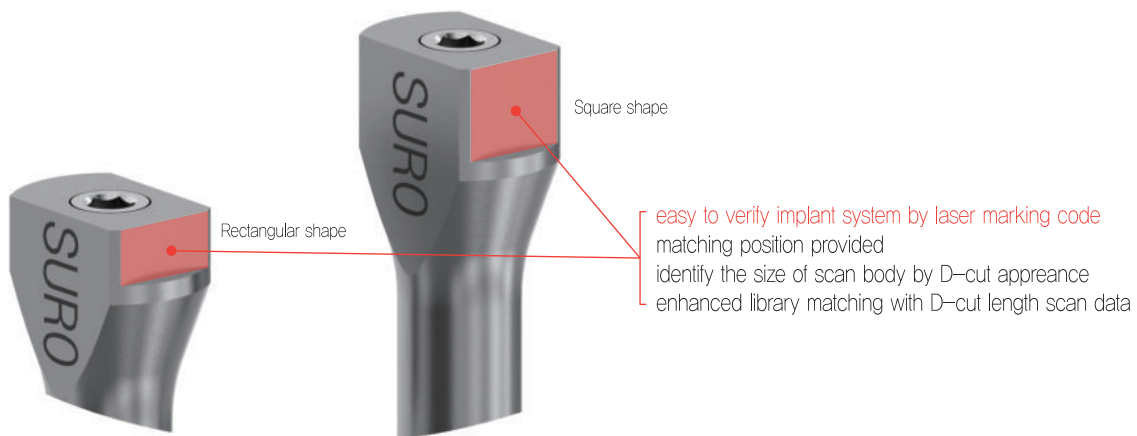
Design concept

TiN Screw

- Surface treatment of TiN-coating
- increased of abrasion resistance
 - improved fixation of screw

Screw made of Titanium Grade ELI
Two Screw – TiN coating screw & Titanium screw
TiN coating screw – Final Setting

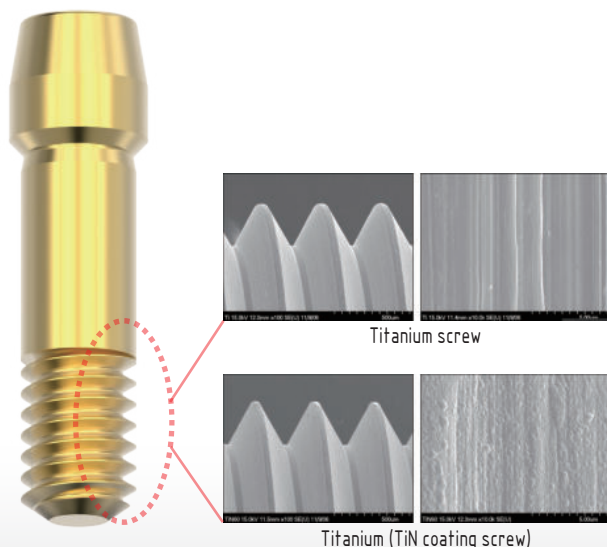
Titanium screw

Titanium (TiN coating screw)

Ti Base

Screw channel path can be changed for different angle
(ball type angled driver necessary)

LS Multibase for Anterior.
Customized Zirconia abutment.
Grooves – Good adhesive strength.
No need for sandblasting.
Tin-coating – High esthetic.

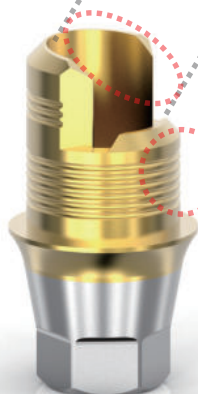Geo Multibase & Screw made of Titanium Grade ELI
Two Screw – TiN coating screw & Titanium screw

Design concept

Ti Base

LL Multibase height is editable on Geo Library.


Tin-coating - High esthetic.
LL Multibase for Posterior.

Grooves - Good adhesive strength.
No need for sandblasting.

Ideal design for making emergence profile.

NEW Item

Gingival height - 2.0mm
height - 6.0mm

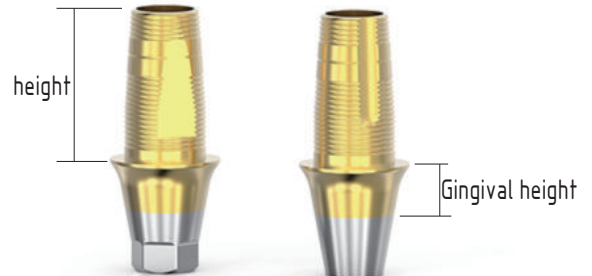Connection

Gingival height-1.3mm --> height-5.5mm

Gingival height-2.0mm --> height-6.0mm

Two different gingival height is available.

Zirconia - Screw retained restorations
Customized Zirconia abutment
Geo Multibase & Screw made of Titanium Grade ELI
Two Screw - TiN coating screw & Titanium screw

Geo Screw Driver

Angled Driver

Square Type / Ratchet Type
 Use with LS Multibase Abutment
 Maximum Angle 25°
 Recommended torque allowance : below 30Ncm.(reference manufacture's guide)
 Caution: unstable connection may bend or destroy the driver
 Recommended usage (1.2hex/1.27hex)-10 times
 (Universal Type)- 20 times
 Connect with Square type Torque wrench
 May use Round Type Torque wrench
 by connecting HDEC - Short (Clinic)Handle

Connect Straight & Angle driver
with the Torque Wrench

Straight Driver

Square Type / Ratchet Type
 Connect by manufacture's recommended torque allowance
 Damage on the tip or usage of long period of time need to exchange driver
 Connect with Sqaure type Torque wrench
 May use Round Type Torque wrench
 by connecting HDEC - Short (Clinic)Handle

Universal - Round hexagon

Square Type / Ratchet Type
 Connect by manufacture's recommended torque allowance
 Damage on the tip or usage of long period of time need to exchange driver
 Connect with Square type Torque wrench
 May use Round Type Torque wrench
 by connecting HDEC - Short (Clinic)Handle

Angled Driver

1.25mm Hex

			
ASDZIL	ASDZIR	ALDZIL	ALDZIR
Square -20mm	Square -13mm	Ratchet -20mm	Ratchet -13mm

CCM Casting Abutment

CCM Casting Abutment

Asia's first developed product in 2009
 Six years of accumulated clinical data
 Excellent biocompatibility
 Incl. screw
 Optimum platform design for emergence profile
 Could make all type of implant restorations

Description

KISPLANT SYSTEM

*set (CCM Abt + screw)

 Single Bridge	Platform Ø Connection	Standard  Ø3.5/4.0/4.5		Driver (Ncm)  1.25Hex
Single	GH : 1mm		DER-H2	25-30Ncm
	GH : 2mm		DER20-H2	
	GH : 3mm		DER30-H2	
Bridge	GH : 1mm		DER-N2	
	GH : 2mm		DER20-N2	
	GH : 3mm		DER30-N2	

Geo


Multi-screw Abutment

- CAD/CAM restoration / Geo Multi-screw Library
Implant Bridge for Zirconia, Titanium and Cobalt Chrome
Zirconia Hybrid
'All On 4" / Hybrid Bar / Bar Overdenture
- Conventional & Digital
- Overdenture, Screw-Retained Restorations.


Multi-screw Abutment & LB-Cylinder

● Multi-screw Abutment

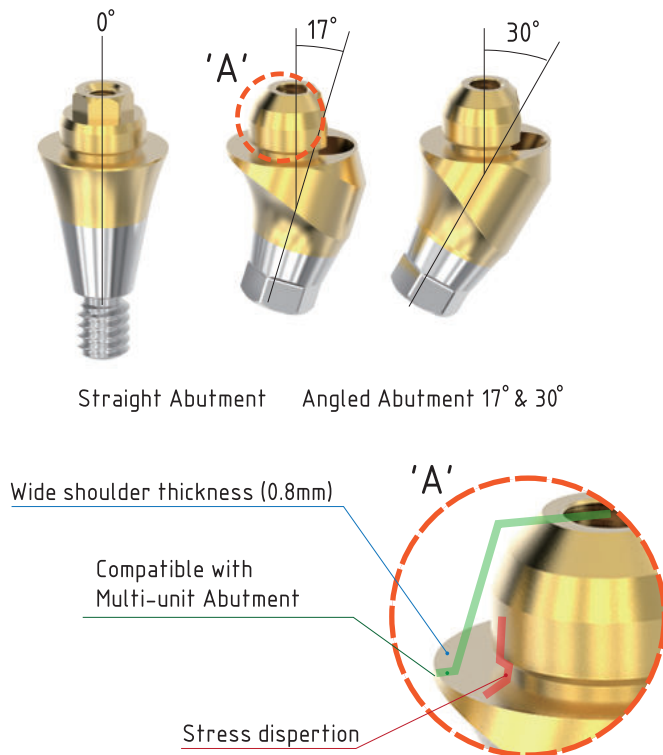● LB-Cylinder

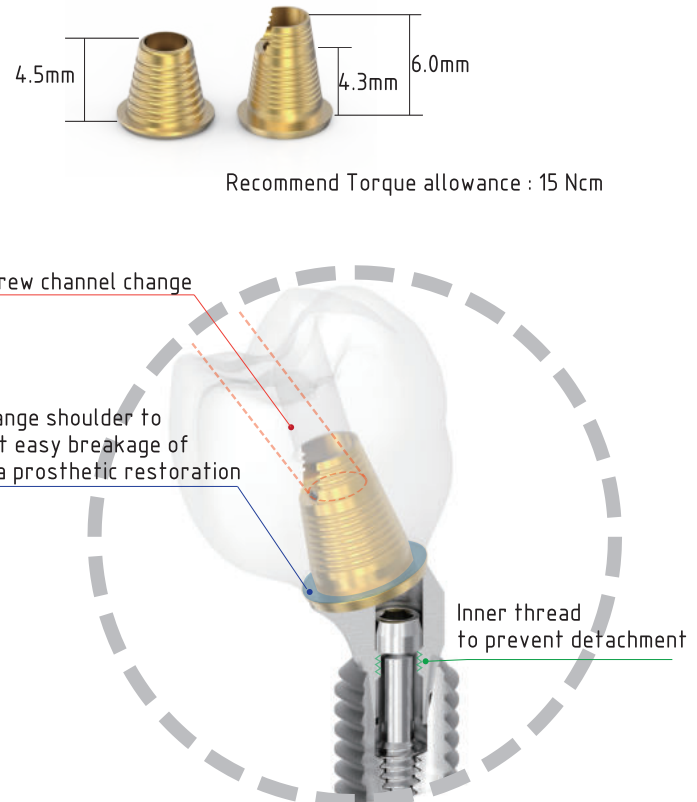Multi-screw prosthetic components of CAD/CAM restoration

● Scanbody

T-Scanbody

Can be used for both scanbody & transfer impression coping
Necessary for angled analog
Necessary for Cylinder screw hole direction

H-Scanbody

Healing Cap is also functional for scanbody with scan spray

SC2-Scanbody


Model scanbody & Intraoral scanbody

● 3D Analog

DT - Series KISPLANT MULTI - UNIT


Multi-screw Straight Abutment

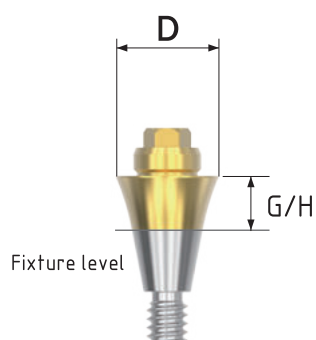Fix with special Straight Guide
Recommended torque allowance : 30Nm (Driver code: G-DV)
packing unit : Multi-screw Straight abutment
+ Straight Guide

G/H		1.5mm	2.5mm	3.5mm	4.5mm	5.5mm
D 4.8						
	Code	G-S15DER-S	G-S25DER-S	G-S35DER-S	G-S45DER-S	G-S55DER-S

Multi-screw Angled Abutment

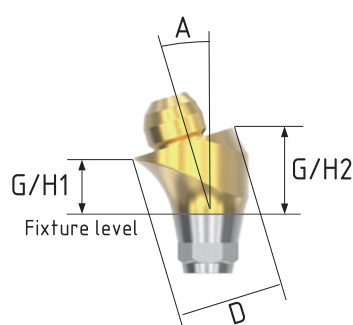Exclusive Angled Abutment Screw
Fix with 1.27 Hex Driver
Implant fixture angle covers up to 108°
Recommended torque allowance : 30Ncm
packing unit : Multi-screw Angled Abutment
+ Angled Guide + screw

Angled Abutment Guide

G/H2		3.0mm	4.0mm	5.0mm	4.0mm	5.0mm	6.0mm	
G/H1		1.5mm	2.5mm	3.5mm	1.5mm	2.5mm	3.5mm	
A		17°	17°	17°	30°	30°	30°	
D4.8								
	Code	G-A1715DER-S	G-A1725DER-S	G-A1735DER-S	G-A3015DER-S	G-A3025DER-S	G-A3035DER-S	

Multi-screw prosthetic components

Product	Healing cap & H-Scanbody	T-Scanbody & Transfer Impression Coping	SC2-Scanbody	Pick-up Impression Coping	3D Analog (Straight)	3D Analog (Angled 17°)	3D Analog (Angled 30°)
Image							
Code	HC-G4840D	IP-G48TD	SC2-G48D	IP-G48PD	3D-G48S	3D-G4817A	3D-G4830A

Product	Temporary Cylinder	LB Cylinder (Short)	LB Cylinder (Long)	CCM Cylinder	Cylinder Screw
Image					 1.25 Hex Driver 15 Ncm
Code	TA-G48D	LB2-G4840D	LB2-G4860D	G48D	MUTDC

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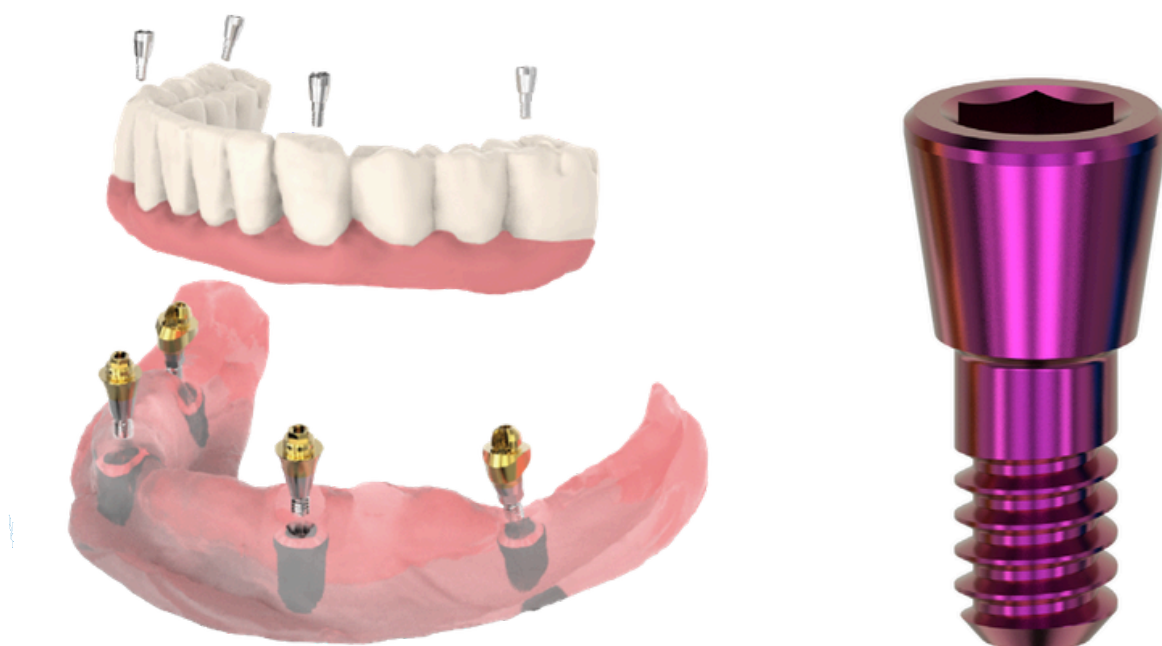
5mm

6mm



Kerator has the world's most powerful competitiveness in overdenture attachment system.



Why?

KERATOR was chosen over the other overdenture attachment, for the comfortable retention power of caps and excellent mobility in the denture metal housing, in the world market.

Clinical Pictures

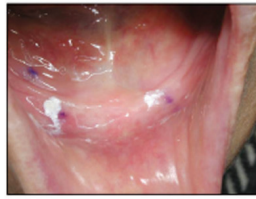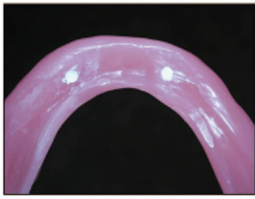
Case A (Mandible)

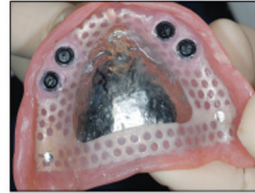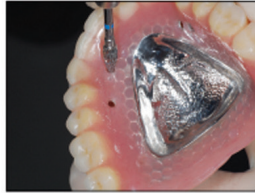
Case B (Maxillary)


Case C (Root Post Type)

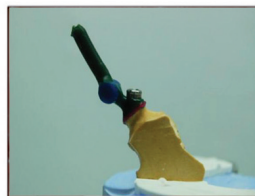
Case D (Root Cast Type)

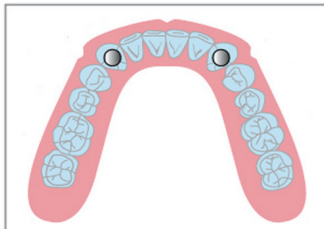
Implant Mandible
Overdenture A

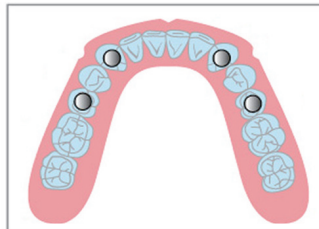Implant Mandible
Overdenture B


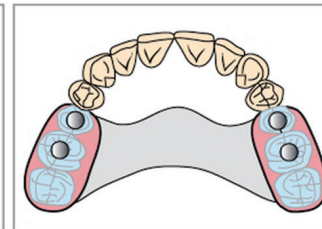
Implant Partial Denture

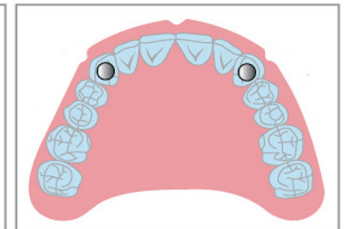Implant Maxilla Overdenture

KERATOR Instruments & Components

KERATOR Magic Tool & Tip, Housing

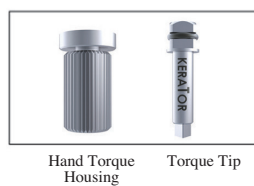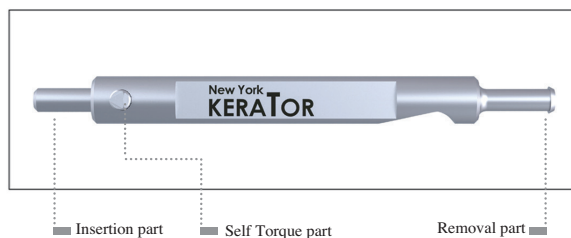Hand Torque
Housing


Torque Tip

KERATOR
Coping

KERATOR
Analog

KERATOR Male Cap Set

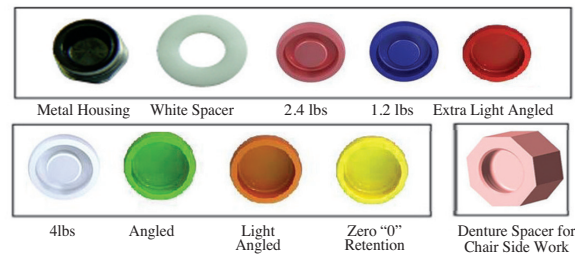Metal Housing

White Spacer

2.4 lbs

1.2 lbs

Extra Light Angled

4lbs


Angled

Light
Angled

Zero "0"
Retention

Denture Spacer for
Chair Side Work

Replacing Male Inserts

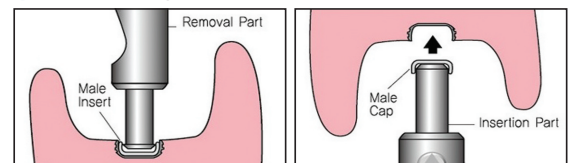
Removal Part

Male Insert


Male Cap

Insertion Part


Fabrication Instructions

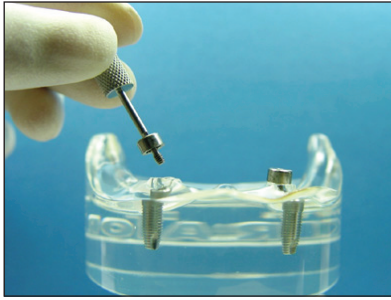
1. Remove healing abutment using hex driver. (0.05inch, 1.27mm)

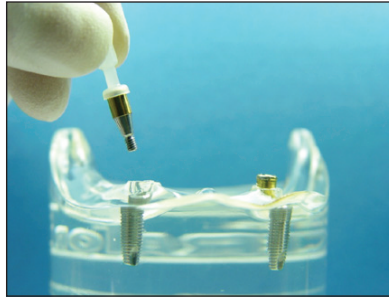2. Engage abutment into the implant using carrier.

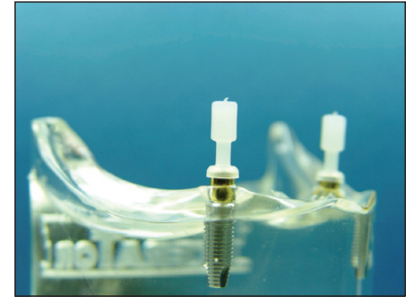
3. Engage straight abutment.
(*check the implant alignment)

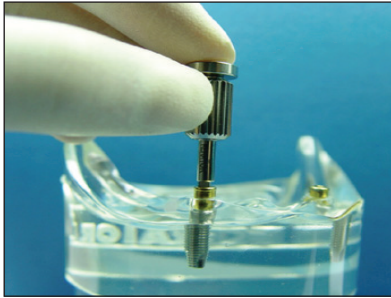4. Tighten the abutment by Hand torque housing & tip


5. Use the Torque Wrench 30Ncm

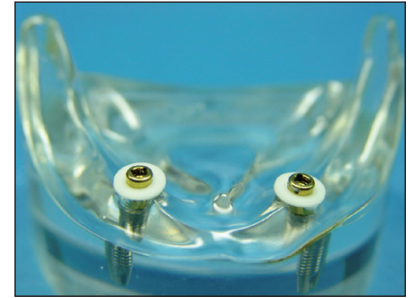6. Place a White Block-out Spacer over the head of each KERATOR abutment.

7. Insert Metal Housing with black cap.

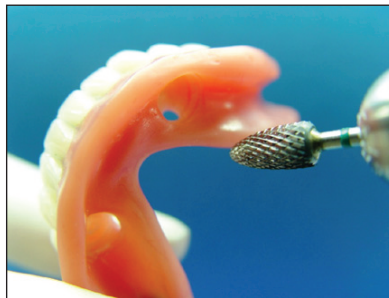8. Prepare a recess in the denture.
Make a vent hole on the lingual aspect to express resin.

9. Mix a resin and place a small amount in the recess of the denture.


10. Insert the denture into position.


11. After the resin has cured, remove the denture and discard the white spacer.

12. Remove the Black cap from the Metal Housing using KERATOR Magic Tool.

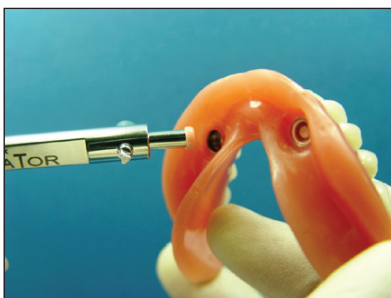13. Push a final cap using KERATOR Magic Tool.

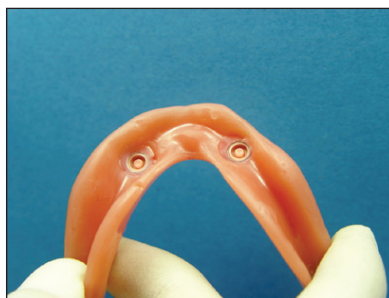
14. Finished figure.


15. Have the patient insert.

Surgical KIT-F

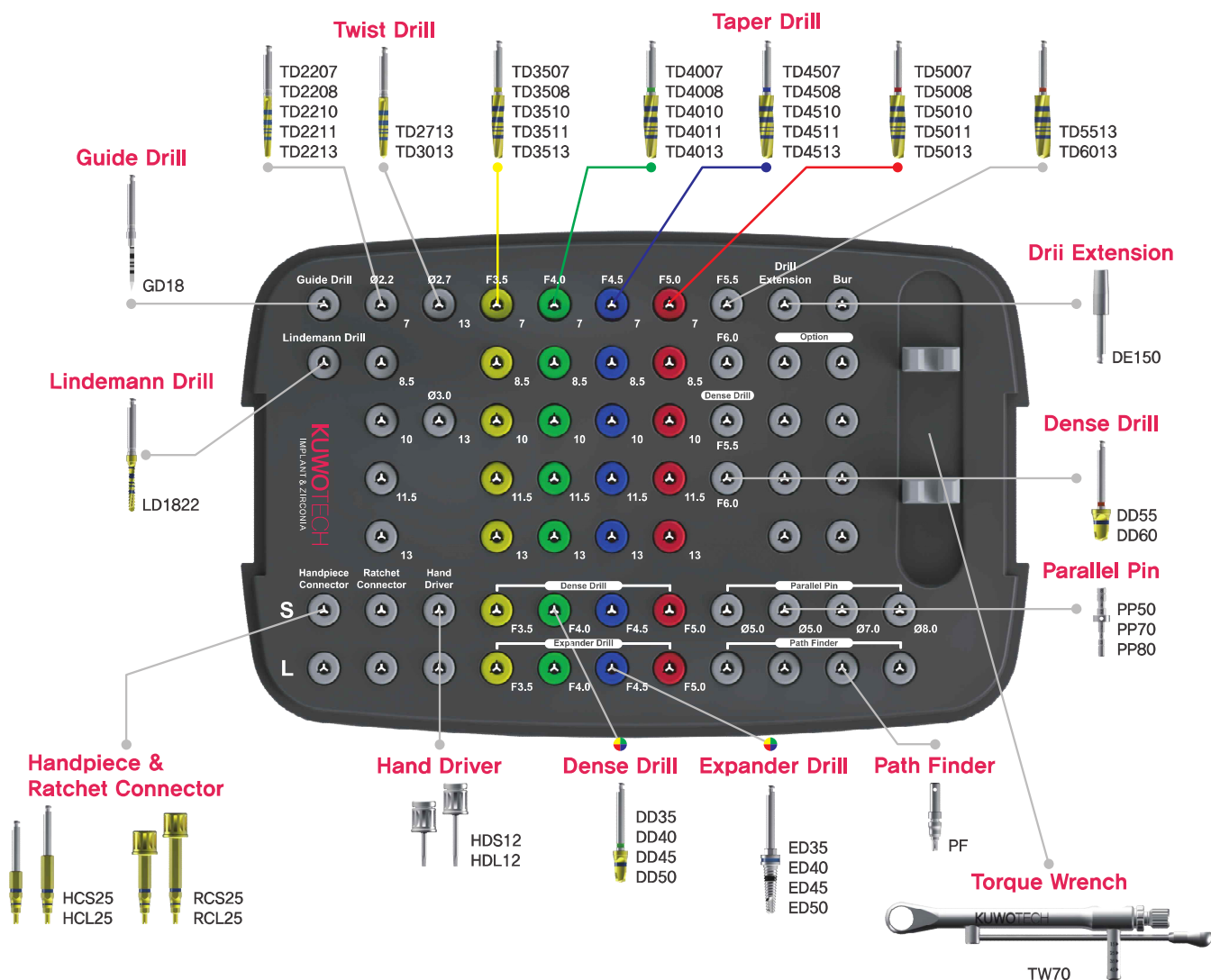Surgical KIT – Stopper KIT (Ref.C STK)

products from KUWOTECH

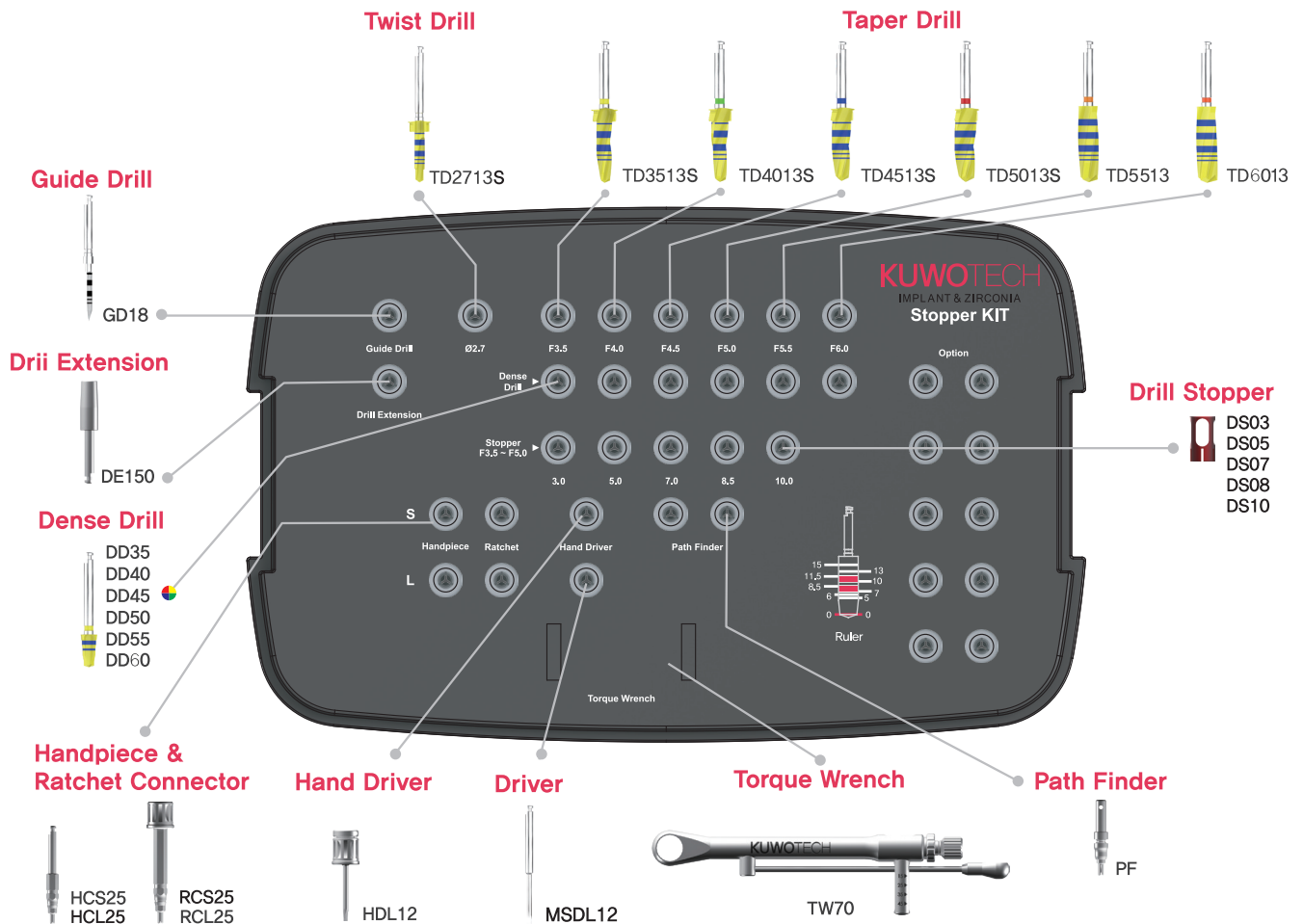Surgical KIT - Drills & Instrument

Guide Drill

· Used to pierce the cortical bone initially

Length	13.0
	
Ø1.8	
Ref.C	GD18

Twist Drill


Length	7.0	8.5	10.0	11.5	13.0
					
Ø2.2					
Ref.C	TD2207	TD2208	TD2210	TD2211	TD2213

Length	13.0
	
Ø2.7	
Ref.C	TD2713


Length	13.0
	
Ø3.0	
Ref.C	TD3013

Taper Drill

Length	7.0	8.5	10.0	11.5	13.0
					
Ø3.5					
Ref.C	TD3507	TD3508	TD3510	TD3511	TD3513

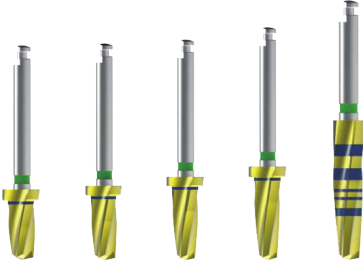Length	7.0	8.5	10.0	11.5	13.0
					
Ø4.0					
Ref.C	TD4007	TD4008	TD4010	TD4011	TD4013

Length	13.0
	
Ø5.5	
Ref.C	TD5513

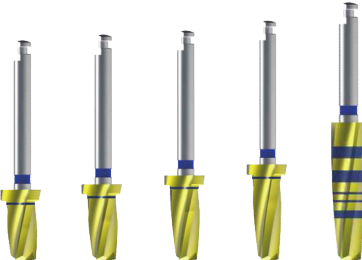Length	7.0	8.5	10.0	11.5	13.0
					
Ø4.5					
Ref.C	TD4507	TD4508	TD4510	TD4511	TD4513

Length	7.0	8.5	10.0	11.5	13.0
					
Ø5.0					
Ref.C	TD5007	TD5008	TD5010	TD5011	TD5013

Length	13.0
	
Ø6.0	
Ref.C	TD6013

Dense Drill

- Used to remove and shape cortical bone to control initial stability in dense bone (Bone type D1/D2)

Diameter	3.5	4.0	4.5	5.0	5.5	6.0
						
Ref,C	DD35	DD40	DD45	DD50	DD55	DD60

Expander Drill

- Useful in very narrow bone

Diameter	3.5	4.0	4.5	5.0
				
Ref,C	DD35	DD40	DD45	DD50

Lindeman Drill

- Can correct the path during drilling

Diameter	1.8 – 2.2
	
Ref,C	LD1822

Drill Extention


- Can correct the path during drilling

Ref.C	DE150
-------	-------

Handpiece Connector


- Delivers torque for the placement of a fixture with a handpiece
- Easy and secure pick up and delivery

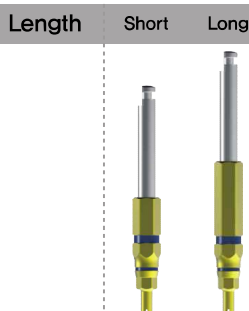Ref.C	HCS25	HCL25
-------	-------	-------

Ratchet Connector

- Delivers torque for the placement of removal of a fixture with a ratchet wrench
- Easy and secure pick up and delivery

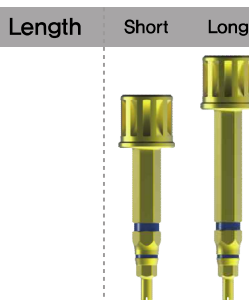Ref.C	RCS25	RCL25
-------	-------	-------

Hand Driver

- Used for all cover screws, all abutment screws and all healing abutment

Ref.C	HDS12	HDL12
-------	-------	-------

Path Finder

- After placing a fixture, a path finder can be connected to guide parallel for the next implant
- Gingival depth can be measured with the grooves especially for flapless surgeries


Ref.C	PF
-------	----

Parallel Pin

- Delivers torque for the placement or removal of a fixture with a ratchet wrench
- Easy and secure pick up and delivery
- Use for checking the direction and location for bone preparation
- Predicts the diameter of an abutment to be secured

Diameter	5.0	7.0	8.0
			
Ref.C	PP50	PP70	PP80

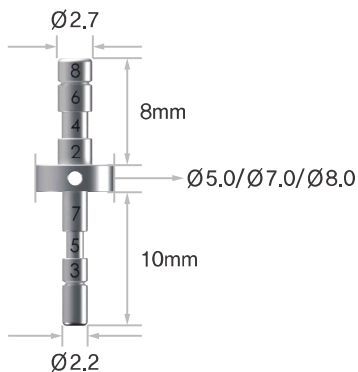Torque Wrench


- Torque wrench has torque option from 15Ncm to 45Ncm and is used for the placement of an implant and final tightening of the abutment screw

Ref.C	TW70
-------	------


Surgical Kit-F Drilling Sequence

BT3510

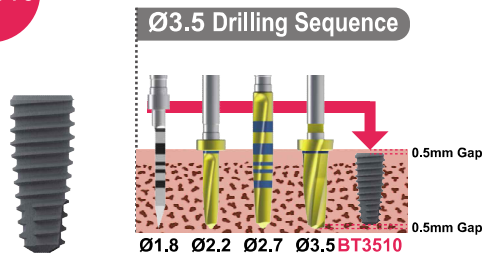
10.0mm is the fixture length and the Taper drill are 0.5mm longer than the fixture, total drill depth is 10.5mm

BT4010

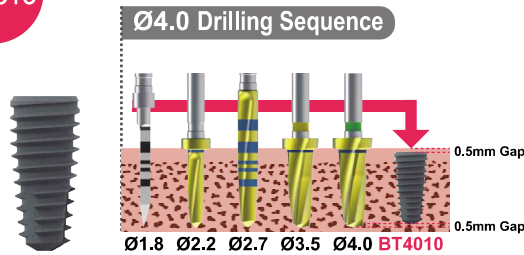10.0mm is the fixture length and the Taper drill are 0.5mm longer than the fixture, total drill depth is 10.5mm

BT4510

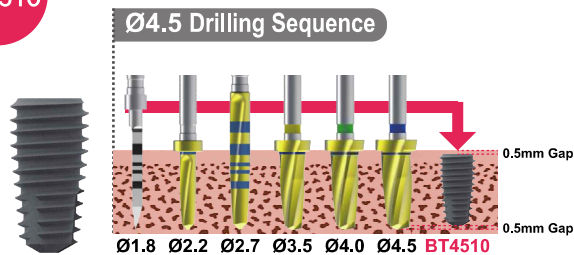
10.0mm is the fixture length and the Taper drill are 0.5mm longer than the fixture, total drill depth is 10.5mm

BT5010

Ø5.0 Drilling Sequence

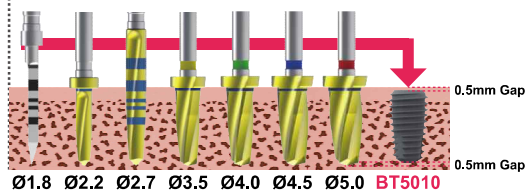
10,0mm is the fixture length and the Taper drill are 0,5mm longer than the fixture, total drill depth is 10,5mm

BT5510

Ø5.5 Drilling Sequence

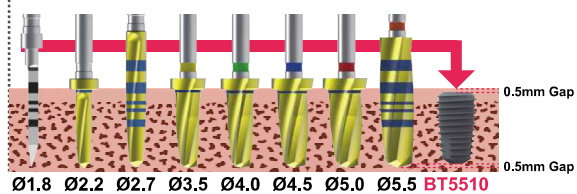
10,0mm is the fixture length and the Taper drill are 0,5mm longer than the fixture, total drill depth is 10,5mm

BT6010


Ø6.0 Drilling Sequence

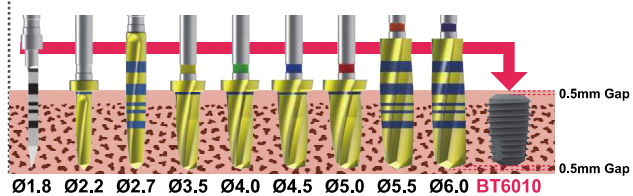
10,0mm is the fixture length and the Taper drill are 0,5mm longer than the fixture, total drill depth is 10,5mm

Taper Drill

- Each drill has depth marking lines from 7.0mm to 15.0mm

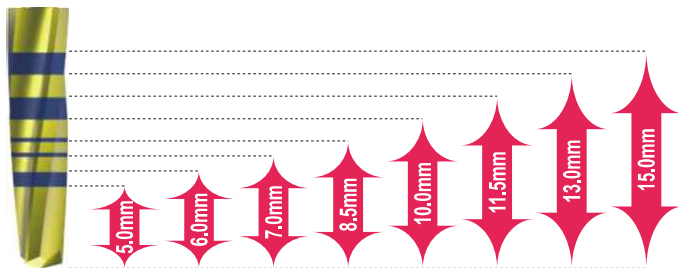
* Taper Drill markings are 0,5mm longer than the fixture so fixtures will automatically be placed 0,5mm subcrestally if the drilling protocol is followed


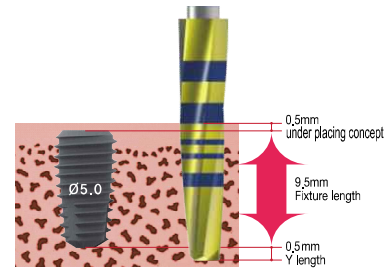
Dense Drill

- To control initial stability in Dense bone (type I & II) use the Dense drill to remove and shape the cortical bone.

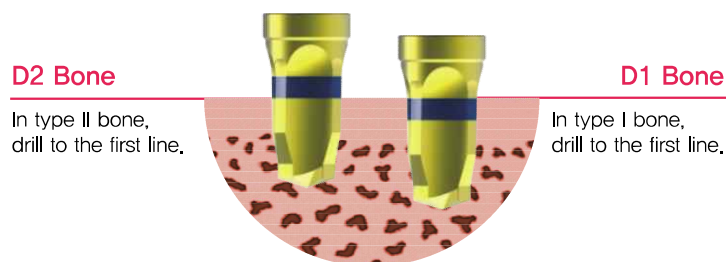
Handpiece & Ratchet Connector

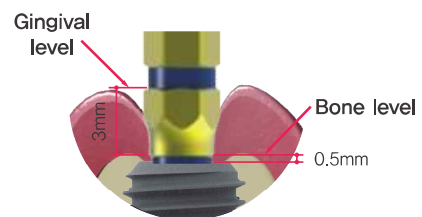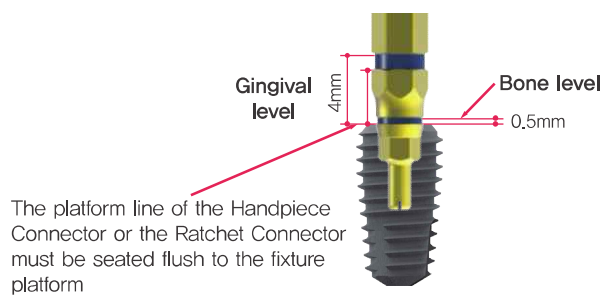
- * Fixture level : Placement should be 0.5mm subcrestal.
- * Gingival level : Line is 3mm above the bone level line and 3.5mm above the platform line

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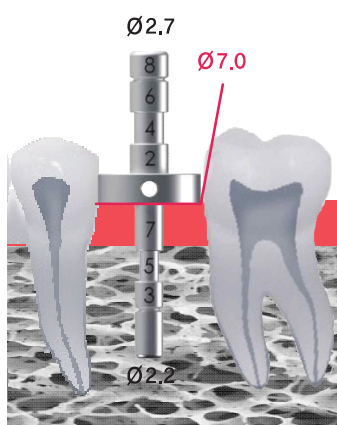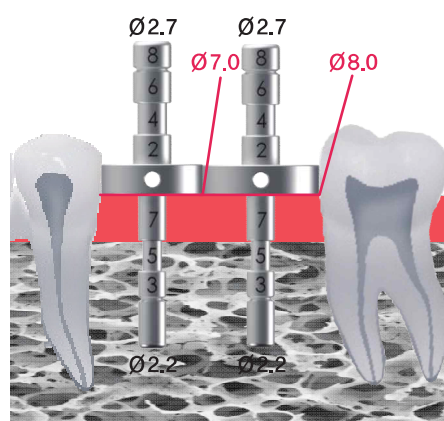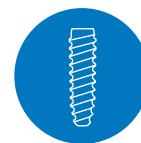
CONNECTION

Platform Switching

Reduces the risk of inflammatory cell infiltration

Greatly contributes to marginal bone preservation and soft tissue attachment

Small but Strong

Suitable for narrow interdental space like maxillary lateral incisors as well as mandibular incisors

Morse Tapered Hex Connection

Provides anti-rotation and allows for convenient indexing of the abutment

Distributes the load to the fixture evenly and helps to reduce the micro movement

Hybrid Double Thread

Allows Easier and quicker implant insertion in various bone condition

Helps to reduce bone stress during insertion and penetrate the bone with ease

3 Blade Cutting Edge


Self-tapping function allows smooth insertion with minimal bone destruction

Rounded implant apex reduces risk of bone perforation


S.L.A Surface Treatment Sandblasted, Large grits, Acid etched

Most researched and proven surface technology

x3.0k

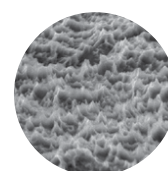x5.0k

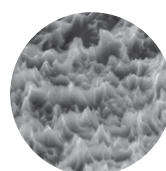Cover Screw

Unit : mm, Scale 1 : 1.2 / mm


Connection	GREEN
-	KMCS14

- Use 1.25 Hex Driver
- Optimum torque : 5-8Ncm
- Screw - M1.4


Healing Abutment

| Diameter Ø3.0

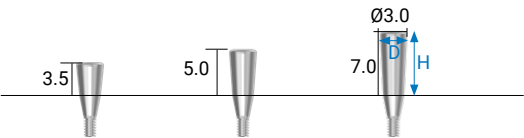Connection	GREEN
H	REF NO.
3.5	KMHA 30 35
5.0	KMHA 30 50
7.0	KMHA 30 70

| Diameter Ø3.5

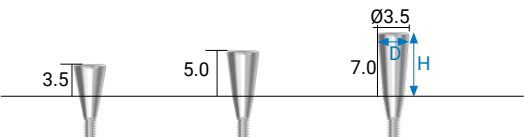
Connection	GREEN
H	REF NO.
3.5	KMHA 35 35
5.0	KMHA 35 50
7.0	KMHA 35 70

- Use 1.25 Hex Driver
- Optimum torque : 8-10Ncm
- Screw - M1.4

Prosthetic Flow Chart Diagrams


Fixture Level

Scale : 1 : 1 / mm

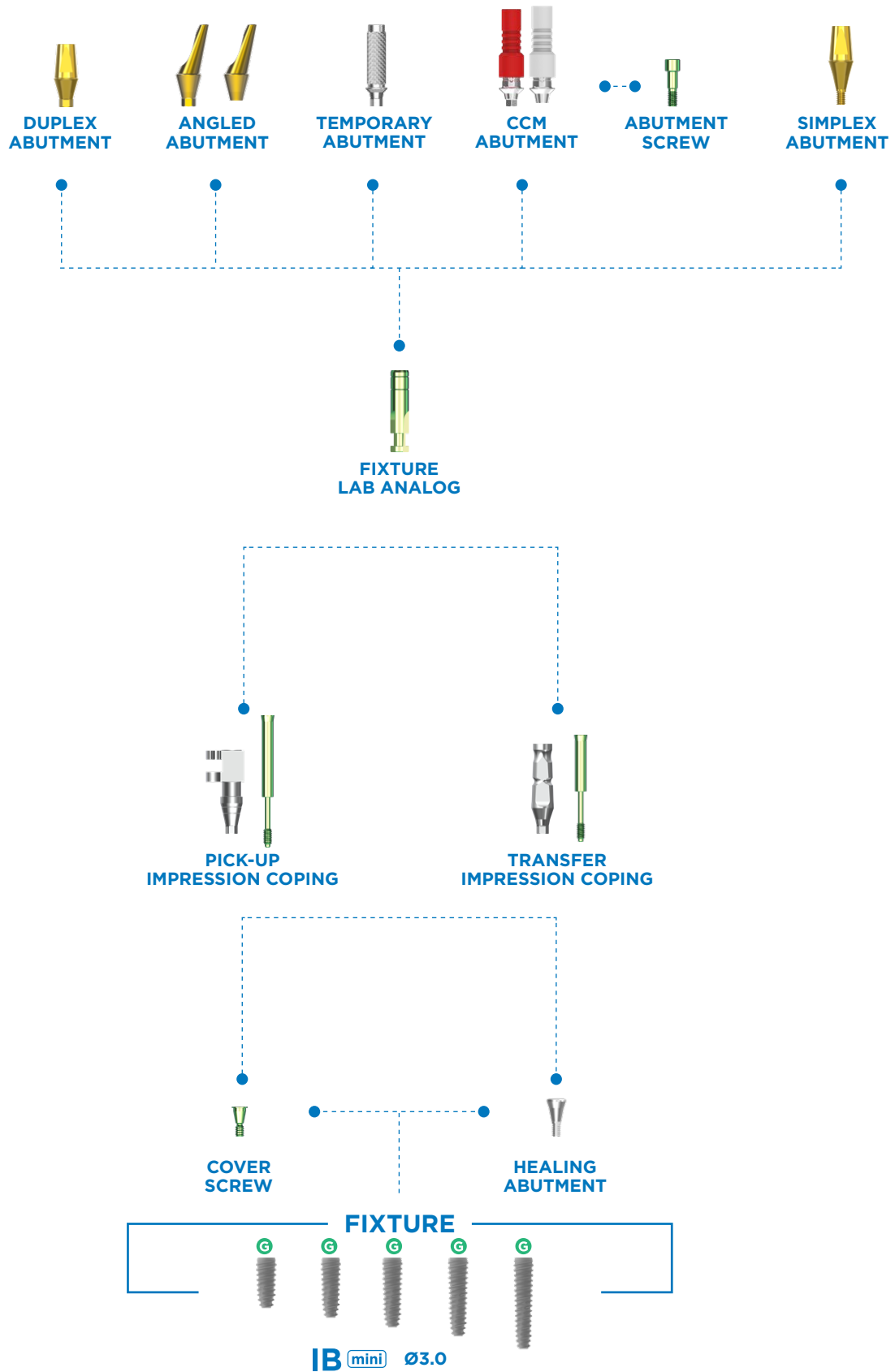
Pick-up Impression Coping

Unit : mm, Scale 1 : 1.2 / mm

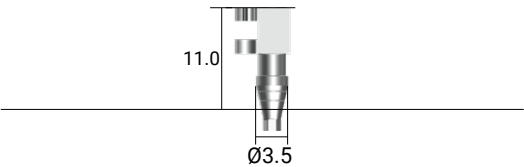
Connection	GREEN
-	KMIP35H


- Pick-up type impression coping with the custom tray
- Hex 1.7 / Screw - M1.4

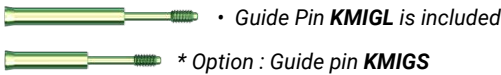Pick-up Impression Coping Guide

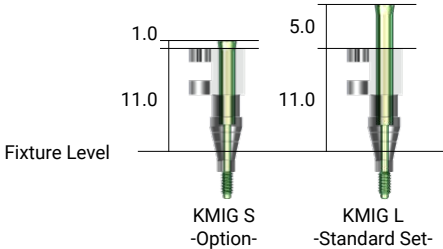
Transfer Impression Coping

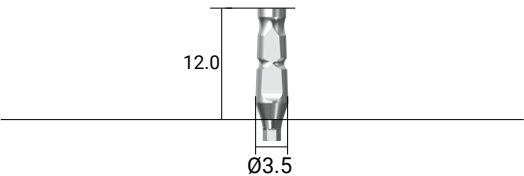
Connection	GREEN
-	KMIT35H

- Transfer type impression coping with the custom tray
- Hex 1.7 / Screw -M1.4

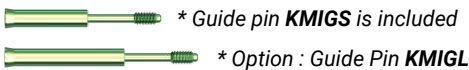Transfer Impression Coping Guide


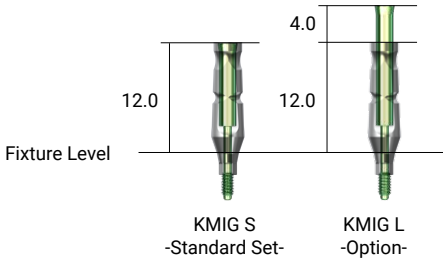
Fixture Lab Analog

Connection	GREEN
-	KMFLA30

- Play a role of fixture on a working model.

Duplex Abutment

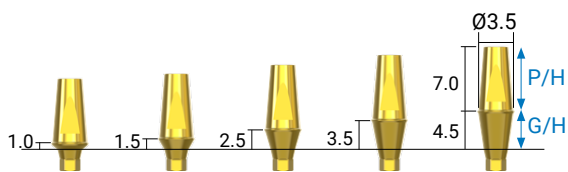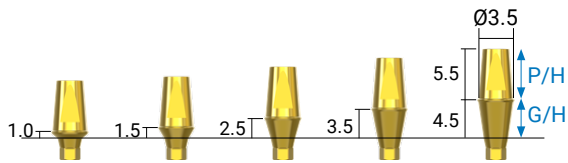Unit : mm, Scale 1 : 1.2 / mm

Connection		GREEN
G/H	P/H	5.5
1.0		KMDA 35 10 55H
1.5		KMDA 35 15 55H
2.5		KMDA 35 25 55H
3.5		KMDA 35 35 55H
4.5		KMDA 35 45 55H

Connection		GREEN
G/H	P/H	7.0
1.0		KMDA 35 10 70H
1.5		KMDA 35 15 70H
2.5		KMDA 35 25 70H
3.5		KMDA 35 35 70H
4.5		KMDA 35 45 70H


- Designed for general cement type prosthesis
- Colored in gold for aesthetic satisfaction
- Use 1.25 Hex driver
- Maximum torque : 25Ncm
- Hex 1.7 / Screw - M1.4


• Abutment Screw **KMAS 14** is included

Angled Abutment

| Angled 15°
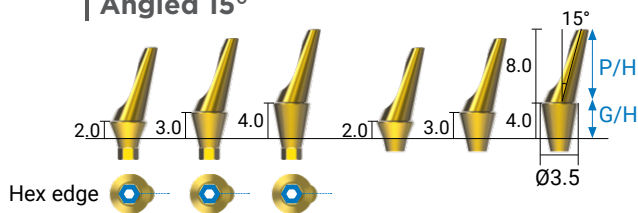
| Angled 25°

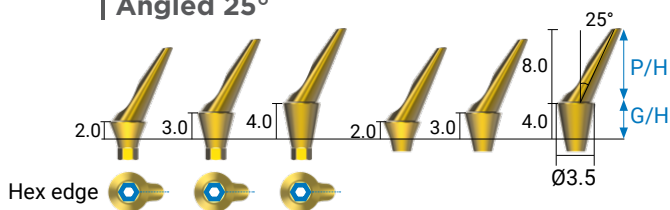Connection		GREEN	
Type		Hex	Non-Hex
G/H	P/H	8.0	
2.0		KMAA 35 2 15H	KMAA 35 2 15N
3.0		KMAA 35 3 15H	KMAA 35 3 15N
4.0		KMAA 35 4 15H	KMAA 35 4 15N

Connection		GREEN	
Type		Hex	Non-Hex
G/H	P/H	8.0	
2.0		KMAA 35 2 25H	KMAA 35 2 25N
3.0		KMAA 35 3 25H	KMAA 35 3 25N
4.0		KMAA 35 4 25H	KMAA 35 4 25N


- It is applied to adjust a path for prosthesis with its 15°, 25° axial angle
- 22° tapered structure for stable connection
- Colored in gold for aesthetic satisfaction
- Use 1.25 Hex driver
- Maximum torque : 25Ncm
- Hex 1.7 / Screw - M1.4

• Abutment Screw **KMAS 14** is included

Temporary Abutment

Unit : mm, Scale 1 : 1.2 / mm

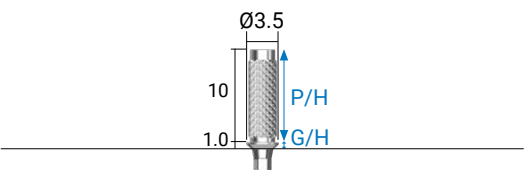
Connection	GREEN
-	KMTA35H

- Use 1.25 Hex Driver
- Maximum torque : 25Ncm
- Hex : 1.7, Screw : M1.4

• Abutment Screw **KMAS 14** is included

CCM Abutment

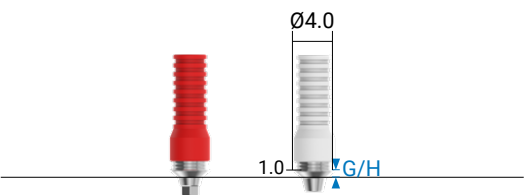Connection	GREEN	
G/H \ Type	Hex	Non-Hex
1.0	KMCC-H	KMCC-N

- Easy manufacturing process due to plastic sleeve
- Improved physical connection strength due to the deep groove on the CCM abutment
- Easy recognition by color on the plastic sleeve
- **Single [red color]** Bridge [white color]
- Use 1.25 Hex Driver
- Maximum torque : 25Ncm
- Hex : 1.7, Screw : M1.4


• Abutment Screw **KMAS 14** is included

Simplex Abutment

Unit : mm, Scale 1 : 1.2 / mm

- Designed for general cement type prosthesis
- One body structure (abutment with screw)
- Colored in gold for aesthetic satisfaction
- Use 1.25 Hex driver
- Maximum torque : 25Ncm
- Screw - M1.4

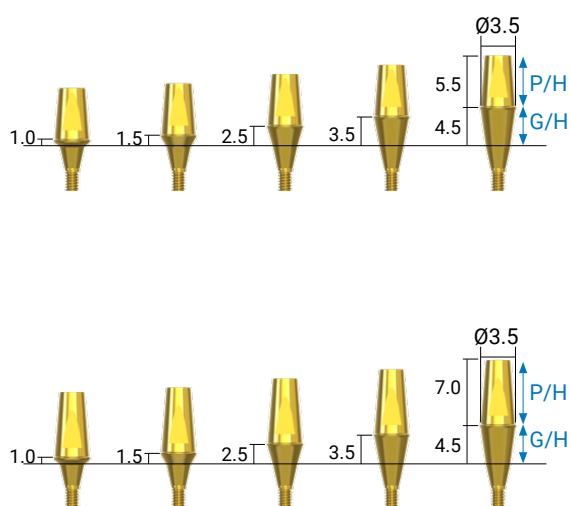Connection		GREEN
G/H	P/H	5.5
1.0		KMSA 35 10 55
1.5		KMSA 35 15 55
2.5		KMSA 35 25 55
3.5		KMSA 35 35 55
4.5		KMSA 35 45 55

Connection		GREEN
G/H	P/H	7.0
1.0		KMSA 35 10 70
1.5		KMSA 35 15 70
2.5		KMSA 35 25 70
3.5		KMSA 35 35 70
4.5		KMSA 35 45 70

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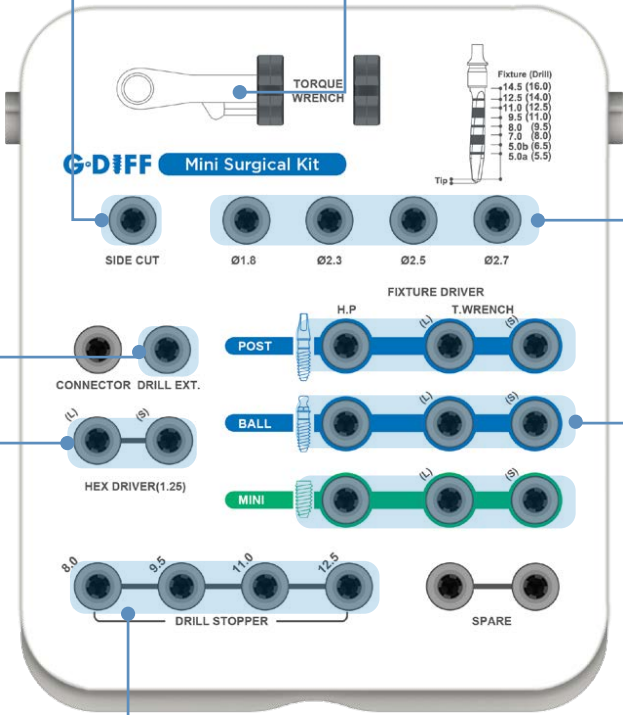
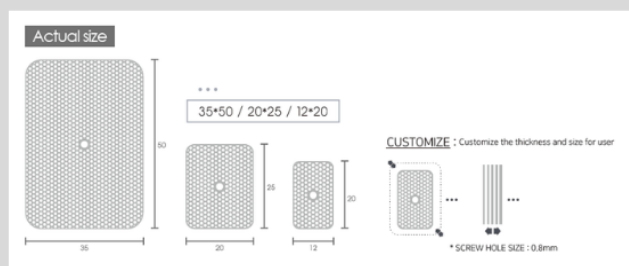
A photograph of two women smiling. The woman in the foreground is a young woman with blonde hair, wearing a light-colored turtleneck sweater. The woman in the background is slightly out of focus, wearing a white jacket with a blue collar. The image is cut off on the right side by a green vertical bar.

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
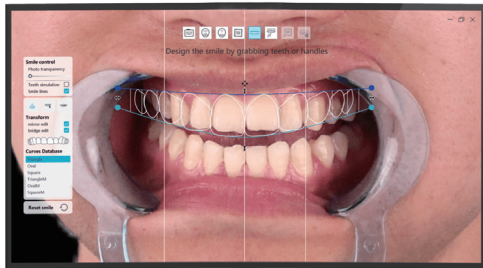
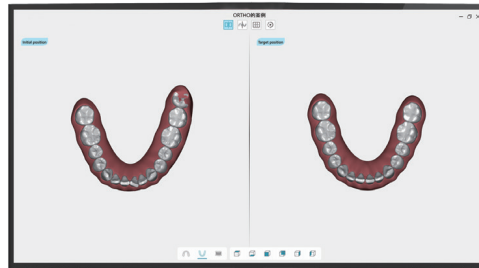
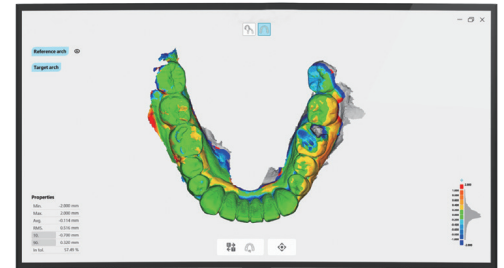
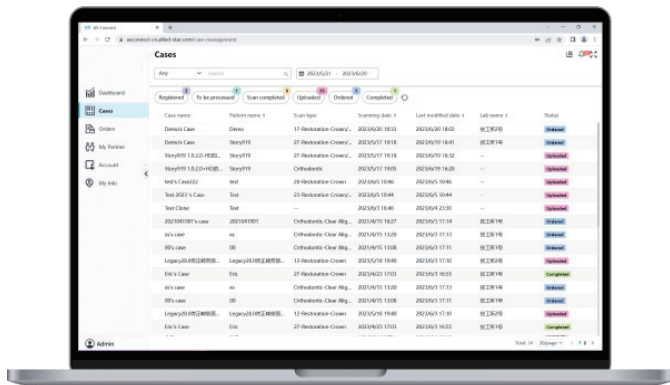
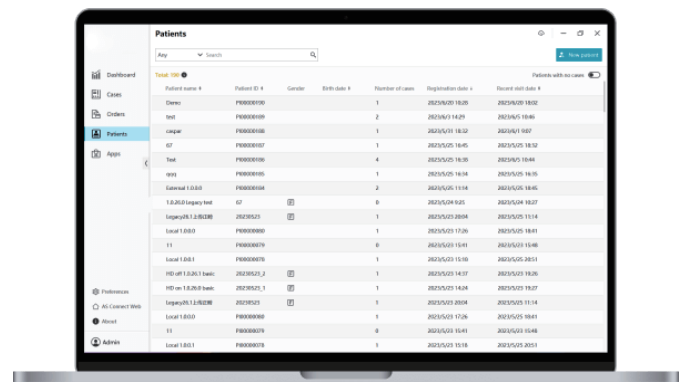
KISPLANT GUIDE KIT

OSSTEM TS3 FULLY GUIDED KIT


- **Fixture Diameter**
Ø3.0, Ø3.5, Ø4.0, Ø4.5, Ø5.0
- **Implant System**
Osstem, Neobiotech, DIO, etc
- **Guided Sleeve**
ALL-GUIDE 5.3 Inner diameter
- **Drill offset**
+9mm

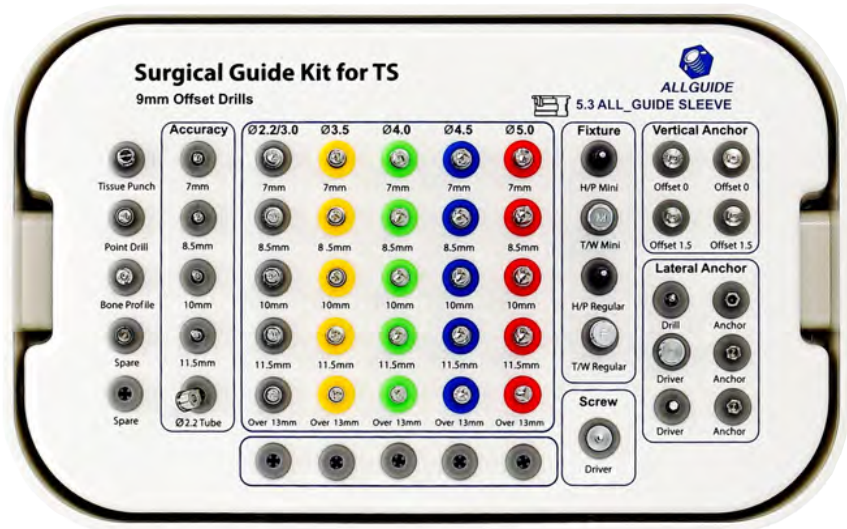Ø5.0 × 7mm + 1.5mm

Sleeve Offset

9mm + 1.5mm

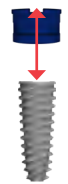※ Please stop
insert fixture
before line

※ Accuracy Option

※ MUST USE Drill

Tissue
Punch

Tube Ø2.2
8.5mm

Point Drill

Ø3.0
7.0mm

Ø3.0
8.5mm

Ø3.5
8.5mm

Ø4.0
8.5mm

Ø4.5
8.5mm

Ø5.0
8.5mm

Fixture
Driver
H.P

Fixture
Driver
Ratchet

※ Depth Control

RPM : 800

RPM : 50 ~ 800

※ Please stop insert
fixture before line

9mm

+1.5mm


KISPLANT GUIDE KIT

KIT COMPONENTS

Point Drill

- A hole on the bone is initially formed on the bone for insertion of the initial drill.


CODE	AGPD
------	------

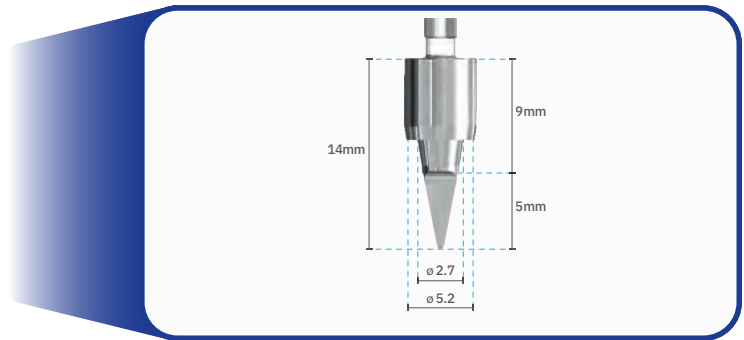Tissue Punch

- Gingiva can be removed stably using a fixed blade inside the punch. (flapless surgery)


CODE	TSTP 3440
------	-----------

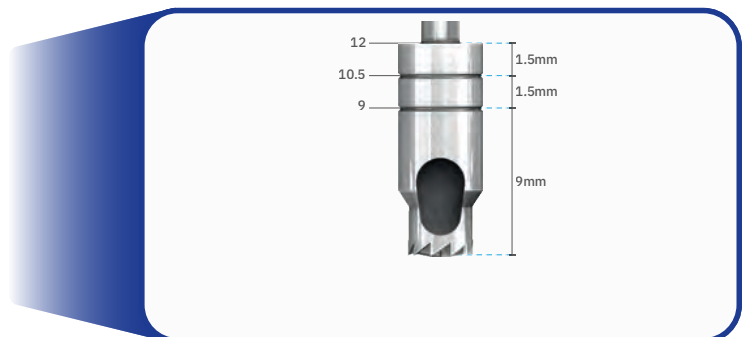Drill tube

- To fix the guide drill and the initial drill without shakage.
- Used with initial drill in sharp ridges.

CODE	AGDT
------	------

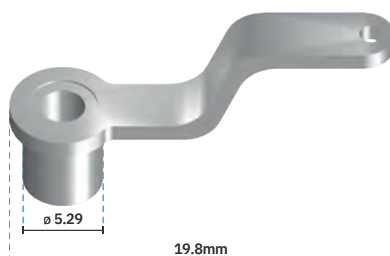KISPLANT GUIDE KIT

KIT COMPONENTS

Initial Drill


- A placement hole(osteotomy site) is formed on the cortical bone.
- 7mm / 8.5mm / 10mm / 11.5mm configuration depending on length.

L \ D	2.2
7	TSID 2207
8.5	TSID 2208
10	TSID 2210
11.5	TSID 2211


Final Drill

- The drill hole is expanded to the final drilling.

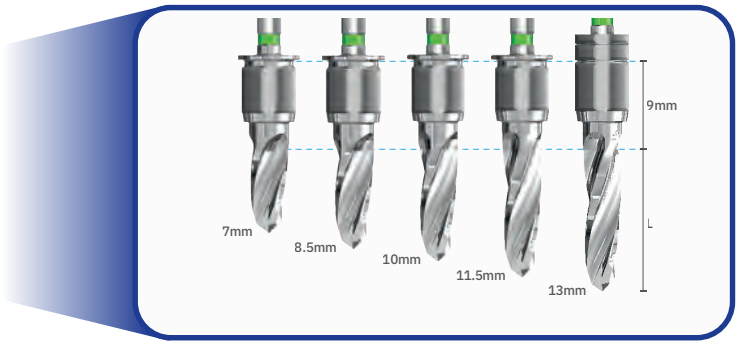- Non-stopper form for 13mm

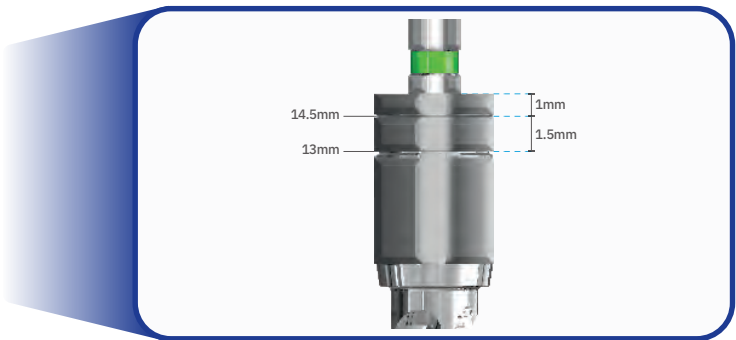L \ D	3.0	3.5	4.0	4.5	5.0
7	TSID 3007	TSID 3507	TSID 4007	TSID 4507	TSID 5007
8.5	TSID 3008	TSID 3508	TSID 4008	TSID 4508	TSID 5008
10	TSID 3010	TSID 3510	TSID 4010	TSID 4510	TSID 5010
11.5	TSID 3011	TSID 3511	TSID 4011	TSID 4511	TSID 5011
13 / 14.5	TSID 3013	TSID 3513	TSID 4013	TSID 4513	TSID 5013


KISPLANT GUIDE KIT

KIT COMPONENTS

Fixture driver (H/P, T/W)

- Implant is placed according to depth of the fixture and direction of the internal hex planned in advance.
- Marked at a 1.5mm Depth, the deeper of the can be adjusted according to the sleeve offset.
- Configured with for Mini Fixture(Yellow code) and Regular Fixture(Green code).

Hand Piece	TSFDHM21
Wrench	TSFDTM21

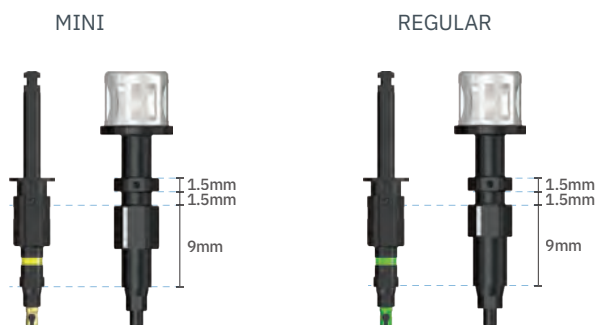Hand Piece	TSFDHR25
Wrench	TSFDTR25

Bone Profile

- The cortical bone disrupting sitting of the abutment is removed

CODE	TSBP
------	------

Hand driver

- Screw driver

CODE	TSHD12
------	--------


KISPLANT GUIDE KIT

KIT COMPONENTS

Anchor drill

- A hole is formed to insert an anchor pin into the bone for the purpose of fixating the surgical guide.


CODE	LAD
------	-----

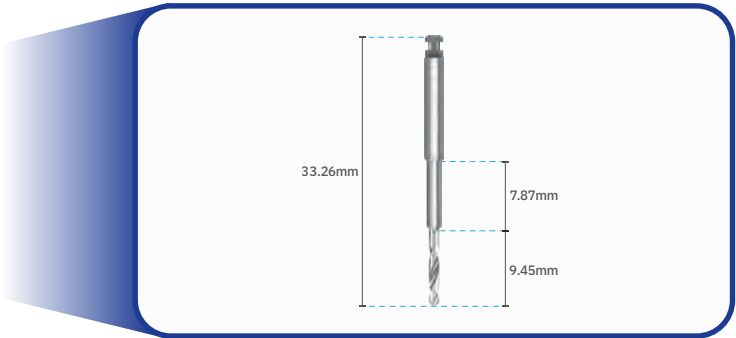Anchor screw

- Fix with surgical guide


CODE	LAC
------	-----

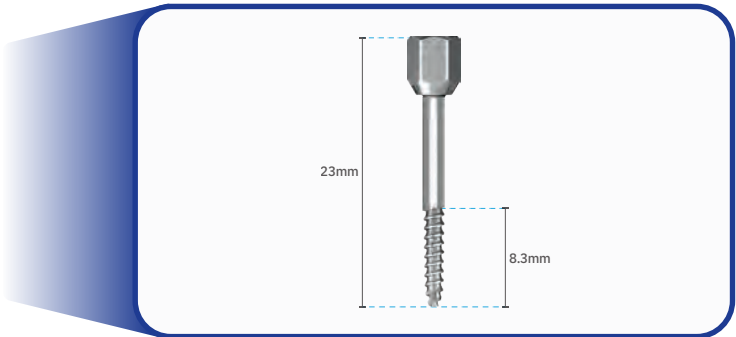Anchor screw driver

- Anchor screw driver for H/P & Wrench

HAND	LAHD
H/P	LAMD

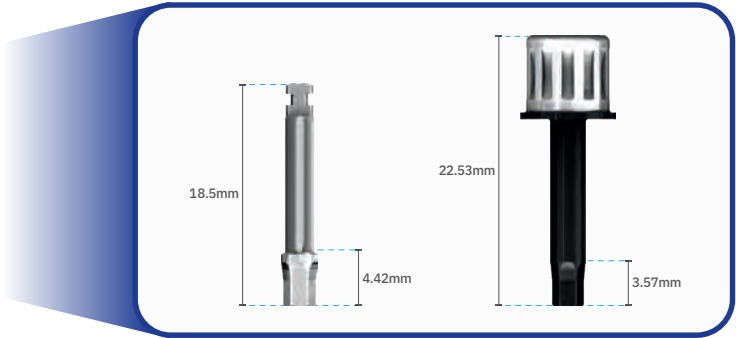
KISPLANT GUIDE KIT

KIT COMPONENTS

Vertical anchor

- For edentulous cases, use it to fix the guide.
- Configured as 0mm, 1.5mm according to offset.

0mm	TSVA 2090S
1.5mm	TSVA 2010L

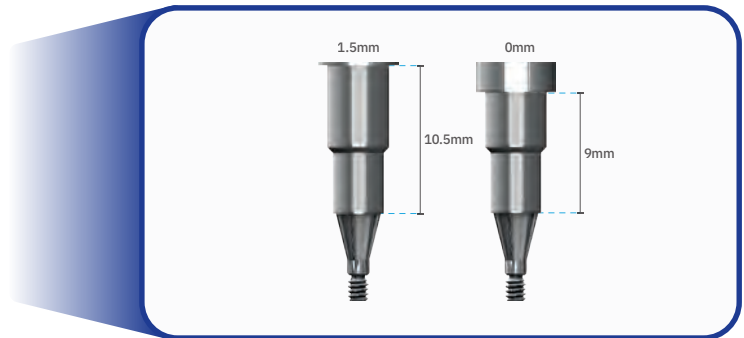Mount connector

- Length of the connector is extended during implant placement.

CODE	TSMD355
------	---------

Torque wrench

- To place the fixture using an implant driver.

CODE	AGTW
------	------